
Elvari Kristallmatte

SKU: ELV-00001
ELV-00009
GTIN: 4170000094326
4262424710057
Datum: September 2023
Version: V8

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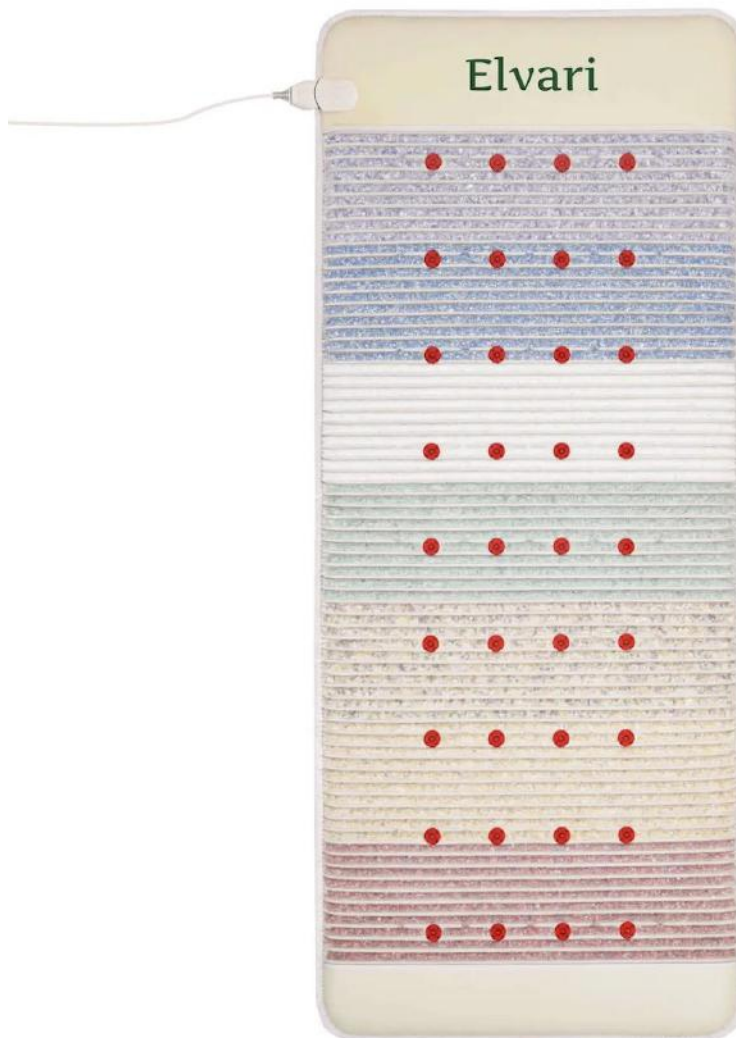English Manual

Our support team will gladly assist you in English as well.
In case you have any questions please contact us via support@elvari.de


Table of contents

Functions	Page 14
First Steps	Page 15
General Notes	Page 16
Programs	Page 17
Recommended programs	Page 19
Questions & Answers	Page 20
Safety Instructions	Page 21
Medical Disclaimer	Page 22
Warranty & Disposal	Page 22
Technical Data	Page 23

Contents of delivery

 <p>Therapy Mat</p>	 <p>Protective Cover</p>  <p>Controller</p>	 <p>Carry Bag</p>  <p>Magnetic Field</p> <p>Hint: Switch on a PEMF program and move the magnetic field tester over the mat. The ball starts to vibrate above the 6 hidden magnetic field coils.</p>
---	--	--

Functions

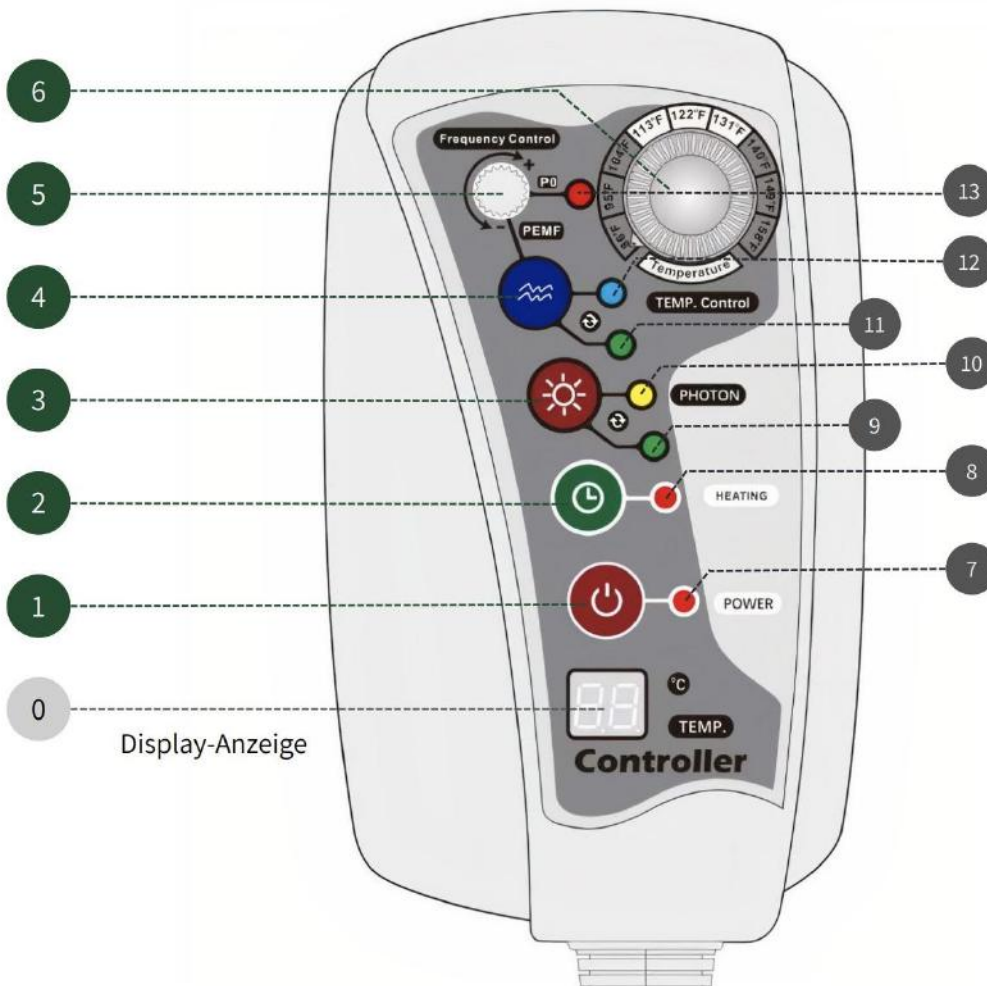Function keys

1. Power (On/Off)
2. Timer button
3. Biophotons button
4. Magnetic field button (PEMF)
5. Magnetic field regulator
6. Temperature regulator

LED indicators

7. Power LED (on = active)
8. Timer LED (on = heating)
9. Biophoton loop LED (lit = P2)
10. Biophoton LED (on = biophotons activated)
11. Green magnetic field LED (on = activated)
12. Blue magnetic field LED (pulses in frequency)
13. Red magnetic field LED (on = manual mode)

First steps

Step I: Unbox the Elvari Therapy Mat

Take the mat out of the carry bag. Place the mat on a flat surface and unfold it to its full length.

Step II: Connect power

Connect the controller to the mat and plug the power plug into an outlet. Turn on the mat by pressing the power button (1) once.

Step III: Set programs

For the beginning we recommend that you use all programs at the same time and start gently at a low level.


Infrared heat: Press the green button (2) twice. The screen now displays "1H" - one hour runtime. Set the temperature to 40 °C via the large scroll wheel (6).

Hint: It takes approximately 10 minutes to warm up all 18 layers.

Biophotons: Press the red button (3) once until the screen shows "P1". The biophotons will now run for 30 min.

Magnetic field (PEMF): Press the blue button (4) a total of 14 times until "PF" is displayed on the screen. This is the program for the Schumann frequency - the heartbeat of the earth.

Important note: The 6 magnetic field coils have different strengths at different locations. This is normal. Read more about this on page 5.

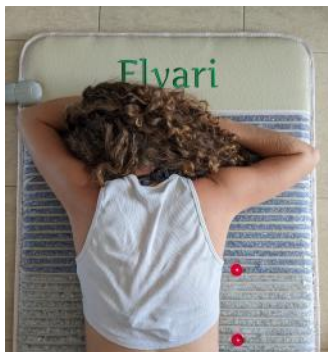Step IV: Relax

Lie down and relax. No cell phone. No laptop. Read a book if you like. Take 30 minutes for yourself.

Our recommendation: Use exactly these settings for 30-60 minutes at least once a day for the first 30 days.

Give your body time to get used to the therapy mat. Then you can start experimenting with the recommended programs (page 6). Listen to what your body tells you.

General Notes

Application Schedule

1. In the first 1-2 weeks, use the recommended programs from "First Steps."
2. Once your body feels ready, switch to a magnetic field program of your choice. Many customers choose based on their symptoms. If you don't have specific symptoms, start with "Detox" (P2) and ideally supplement with Zeolite & Chlorella.
3. After 6-8 weeks, your body will be largely accustomed to the application. Begin experimenting with the intensive programs.

General Notes

1. Drink a full glass of water before and after using the infrared deep heat and/or magnetic field therapy. This is crucial for cellular metabolism.
2. Start your treatment at short and regular intervals. Gradually increase the duration and temperature, focusing on your physical well-being. Health is a marathon, not a sprint!
3. Feedback from our customers indicates that optimal effects can be achieved with 2-4 sessions of 30-60 minutes each week.
4. Try sleeping on/with your device. If it's too hard for you, it's better to use it during the day.
5. Too heavy? Fold your mat and put it in the carry bag. With the handle, it's easier to transport. Alternatively, check out the "Mini Mat" by Elvari.
6. Biophotons and infrared work best on the body when there is as little fabric as possible in between.

Care Tips

1. Protect your device from dirt & stains by using a topper, towel, or even a duvet cover.
2. For cleaning, use a slightly damp cloth. Carefully unplug the device from the socket beforehand.
3. Fold your device only at the designated spots, with the surface facing inward.
4. For infrared sessions above 65°C, we recommend refraining from adding a magnetic field program, as this can impair the magnetic field function.

Program settings

The negative ions and chakra stones work at any time (even without power!) and therefore do not need to be switched on separately. You can easily end programs early by pressing the power button (1).

Temperature settings

You can determine the duration via the green button (2):

1H = 1 hour | 4H = 4 hours | 8H = 8 hours | 0H = 12 hours

You can also use the therapy mat overnight. We recommend a maximum temperature of 40 degrees and that you drink 1/2 liter of water beforehand.

Our recommendation

35 °C - 40 °C: Gentle. Suitable for long periods of use.

40 °C - 50 °C: Standard. Recommended for a duration of 60-120 min.

50 °C - 60 °C: Strong. Suitable for a duration of 10-90 min.

60 °C - 70 °C: Intense. Should be used for a maximum of 10-60 minutes.

i Note on the heating process

After you have switched on the temperature, it can take up to 10 minutes until the mat is completely heated up. This is because the infrared layers are heated first. To speed up the process, you can temporarily turn the temperature up to 70 degrees while the mat heats.

i Note on the thermal protection system

The controller has a safety software that protects you from overheating and dehydration. If the temperature is above 50 °C for more than 3.5 hours, the protection function automatically lowers the temperature to 45 °C and the display shows [].

Biophoton Settings

You can choose between two biophoton programs of different lengths. Program 1 (P1) lasts 30 minutes and lights up yellow. Program 2 (P2) lasts a total of 8 hours and switches on and off every 30 minutes. The LED (9) lights up green during this time.

Hint: The effect of biophotons is only marginally influenced by cover & clothing.

Magnetic field programs

The controller includes a selection of 14 magnetic field programs, all of which have different effects:

P0 Manual program 1-30Hz (7h 20min)	P7 Restful sleep (7h 10min)
P1 Defenses (40min)	P8 Regeneration (7h 6min)
P2 Detox (60min)	P9 General strengthening (60min)
P3 Concentration (60min)	PA Muscle relaxation (60min)
P4 Relaxation (60min)	PC Energy (60min)
P5 Body-Mind-Pain Separation (60min)	PE Awareness (60min)
P6 Connective tissue (60min)	PF Schumann Resonance (76min)

i Notes on the magnetic field LEDs

In the preset modes, the blue magnetic field LED pulses in the current frequency of the magnetic field. These programs also include pauses, during which the blue LED isn't lit. However, you can tell that the program is running properly by the green LED.

You can also adjust the magnetic therapy manually:

Use the program (P0) to enter the manual mode. Now use the scroll wheel (5) to set the exact frequency in Hertz (Hz):

01 Hz Relaxation	09 Hz Body Balance
02 Hz Endorphins	10 Hz Mind-Body Connection
03 Hz Regeneration During Sleep	11 Hz Mindfulness
04 Hz Perception & Peace	13 Hz Vitality
05 Hz Awareness	15 Hz Mind-Body-Pain Separation
06 Hz Feeling Of Security	20 Hz Energization
07 Hz Mental Performance	25 Hz Self-Worth
08 Hz Concentration	30 Hz Dental Health

i Hint: The pulsating magnetic field is strongest directly above the 6 PEMF coils.

To test it, place the magnetic field tester over one of the coils shown on the right. It's normal that the effect is less noticeable in other places.

Many users and alternative practitioners attribute medical effects to the programs, which we can neither confirm nor promise. Please consult your medical specialist if needed.

Recommended programs

In our experience, there are some particularly effective program combinations. If one of them takes too long for you, you can end it prematurely at any time.

<p>Defense Infrared: 35 °C - 45 °C Magnetic field: P1 Biophotons: P1 Application time: 120 minutes</p>	<p>Detox Infrared: 65 °C - 70 °C Magnetic field: P2 Biophotons: P1 Application time: 10 - 20 minutes</p>
<p>Mental peace Infrared: 35 °C - 40 °C Magnetic field: P5 Biophotons: P1 Application time: 30 - 60 minutes</p>	<p>Connective tissue Infrared: 45 °C - 50 °C Magnetic field: P6 Biophotons: P1 Application time: 60 minutes</p>
<p>Bones & joints Infrared: 55 °C - 60 °C Magnetic field: PA Biophotons: P1 Application time: 15 - 90 minutes</p>	<p>Regeneration (night program) Infrared: 35 °C Magnetic field: P8 Biophotons: Off Application time: 6 - 8 hours</p>

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Questions & Answers

> How is the Elvari crystal mat insulated against EMF/E-Smog?

The crystal mat uses a non-EMF heating wire, which has specially reinforced insulation that protects against electromagnetic radiation (EMF). Additionally, a carbon layer is installed, which shields electromagnetic radiation and directs it into the ground.

For those especially sensitive to electricity, we recommend wrapping some aluminum foil around the plug to further insulate it. More information can be found at www.elvari.de/esmog.

> How often and for how long should one use the crystal mat?

We generally recommend a treatment duration of 30-60 minutes per day, every other day. However, one's own bodysense should determine the duration and frequency with which the different therapy forms are used.

Additionally, the mat can also be used for sleeping at night.

> Are there any side effects or contraindications?

In the following cases, we advise against using magnetic field therapy:

- Pregnancy
- If using magnetic/electronic implants (like pacemakers!)
- Epilepsy & seizure-like nerve disorders
- Heart rhythm disorders and heart diseases
- Hyperthyroidism
- Angina Pectoris

All other therapy forms can be set and used independently of the pulsating magnetic field.

Occasionally, there might be an "initial aggravation." This refers to a intensification of the symptoms before the desired improvements occurs. This "initial aggravation effect" is also known from other therapy forms and, by definition, represents the removal of blockages.

Safety instructions

Please read the following safety instructions carefully. Your safety can only be guaranteed if the device is used as intended.

Caution: Risk of electric shock due to water penetration

- Unplug the power cord when the device is not in use.
- Do not pour any liquids over the device.
- Unplug the device before cleaning it.
- Only clean the device with a dry or slightly damp cloth.
- Do not use the device outdoors.
- The power plug and control panel must not come into contact with moisture to avoid the risk of electric shock. Make sure that the product is dry before use.

Caution: Warning of electric shock due to damage of the basic unit or control section

- Regularly check the basic device, controller and accessories for any damage.
- Never use the device with a damaged controller.
- Please do not attempt to repair the device yourself if it is damaged.

Attention: Fire hazard

- Please monitor the temperature carefully while using the device.
- Never use the device folded or crumpled.
- The device can be used with a cotton sheet, but not with thick padding.
- Do not use near fire, heaters or flammable objects.
- Excessive heat can cause skin irritation or burns. Please take care while using the device and lower the temperature if you get too hot.

Attention: General hazard warnings

- Keep the device out of reach of children and let them use it only under supervision.
- Pregnant women and people with pacemakers or magnetic implants must not use pulsed magnetic field therapy.
- If you are affected by a disease associated with adrenal insufficiency, you should not use high temperature settings.
- Limit high temperature settings to a maximum duration of 60 minutes.
- Do not use high temperatures overnight.
- Elderly people and small children should only use the lower temperature settings.

Medical disclaimer

We hereby emphasize that we cannot promise any healing effects with the use of our products, devices or applications and we do not want to do so.

The modes of action mentioned are all based on experience gained during practical application. However, these cannot be guaranteed to the individual case.

For more information, visit www.elvari.de/en/pages/heilmittel-disclaimer

Warranty & Disposal

The legal warranty of 2 years applies, as described in the Consumer Protection Act. If you experience a problem with your device, please contact our support team as the first step.

For more information, visit www.elvari.de/en/pages/30-tage-garantie

Disposal

This symbol on the product or its packaging indicates that the product should not be disposed of as ordinary household waste. Please take your discarded device to a collection point for the recycling of old electrical and electronic equipment.

By properly disposing the therapy mat, you help us to prevent any potential negative impact on the environment. For more information about the disposal and recycling of the product, please contact the relevant municipal waste disposal office, the municipal enterprise or write to our support team.

Declaration of Conformity

The conformity of the product with the legally prescribed standards is guaranteed. The complete declaration of conformity can be requested at the following address:


AMBTS Consulting Ltd, Avgoustou Pikar 14, 3101 Limassol, Cyprus | support@elvari.de

Technical data

Elvari Therapy Mat	190*80*3cm
Temperature range	30-70 °C
Magnetic field strength	150 µT
Weight	16 kg (18 kg shipping weight)
Package dimensions	82*65*15cm
Power consumption	AC 280W 220V 50Hz
Protection type	<input type="checkbox"/> Class 2 Type B Protective insulation

If you have any questions, please contact us at support@elvari.de. We usually reply within 24 hours on weekdays.

Elvari

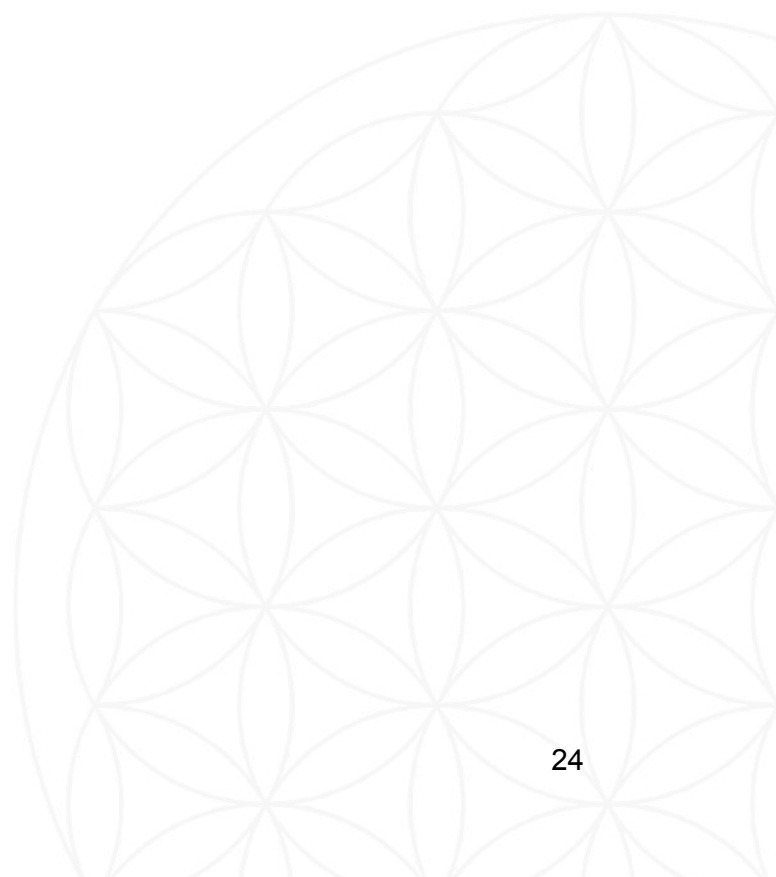Elvari

© Elvari
Eine Marke der AMBTS Consulting Limited
Avogustou Pikar 14, 3101 Limassol, Zypern
support@elvari.de